


VP-3: CARD CAGE ADJUSTMENT (pg 1 of 2)

Time: 10 minutes

Purpose: To ensure customer's ability to insert Assay and System Modules securely

NOTE: To achieve proper adjustment, it may be necessary to remove spacers or adjust the 3 screws mounting the hinge to the Card Cage.

1. Verify that 3 spacers are properly placed under Card Cage:
 - a. Tilt back Card Cage (RR-8.0).
 - b. Position spacers (Figure VP3-1).
 - c. Glue spacers in place.

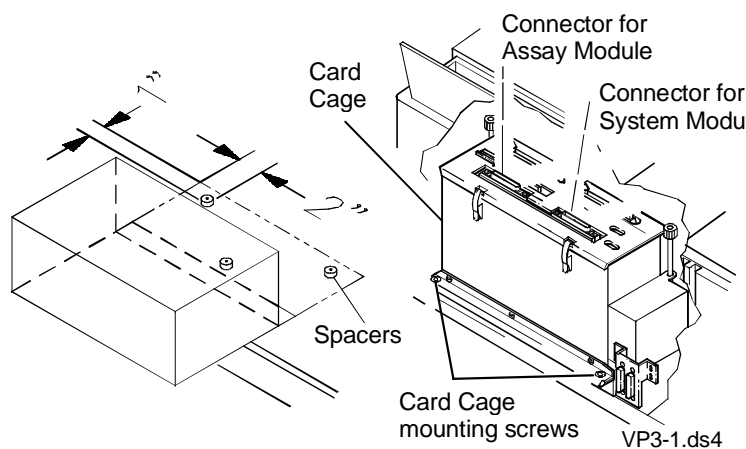

VP-3: Card Cage Adjustment (pg 2 of 2)

2. Loosen Card Cage mounting screws.
3. Adjust Card Cage so that Module connectors on CPU Board align with slots in Enclosure ([Figure VP3-2](#)).
4. Tighten Card Cage mounting screws.
5. Recheck alignment.

