
VP-22: MEIA LAMP ADJUSTMENT (pg 1 of 1)



Time: 2 minutes

Purpose: To align lamp for peak intensity

- 1. Insert MEIA lamp into lamp housing.
- 2. Align either notch with alignment pin.
- 3. Tighten Allen $^{\text{TM}}$ screw.

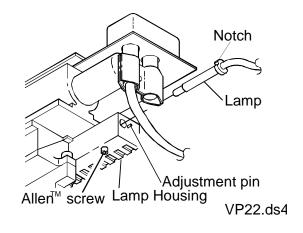

5	Verification Procedures	VP-22: MEIA Lamp Adjustment	5 - 57
	LCx [®] Analyzer Service Manual (Text P/N 47278-102)		

Chapter	Manual	CD Main	Index
Contents	Contents	Screen	