



## RR-8.0: TILT BACK CARD CAGE (pg 1 of 2)

### Removal

1. Remove Enclosure (RR-1.3).
2. Tilt back Card Cage:
  - a. Loosen long captive post screws at front corners of Card Cage.
  - b. Disconnect 4 cables.
  - c. Disconnect Bar Code Cable P8 at Analog Board. (P8 not visible in illustration.)
  - d. Remove Printer cable from cable clamp.
  - e. Tilt Card Cage back. (Refer to [PL-8 illustration](#) in Chapter 3 to see Card Cage in tilted-back position.)


### Replacement

#### CAUTION !

Ensure that no cables are under the captive posts before securing cables.

1. Reposition and secure Card Cage:
  - a. Tilt Card Cage forward.
  - b. Secure Printer Cable in cable clamp.
  - c. Connect Bar Code Cable P8 at Analog Board.
  - d. Connect 4 cables.
  - e. Tighten the long captive post screws to secure the Card Cage to the base.
2. Install Enclosure.

## RR-8.0: Tilt Back Card Cage (pg 2 of 2)

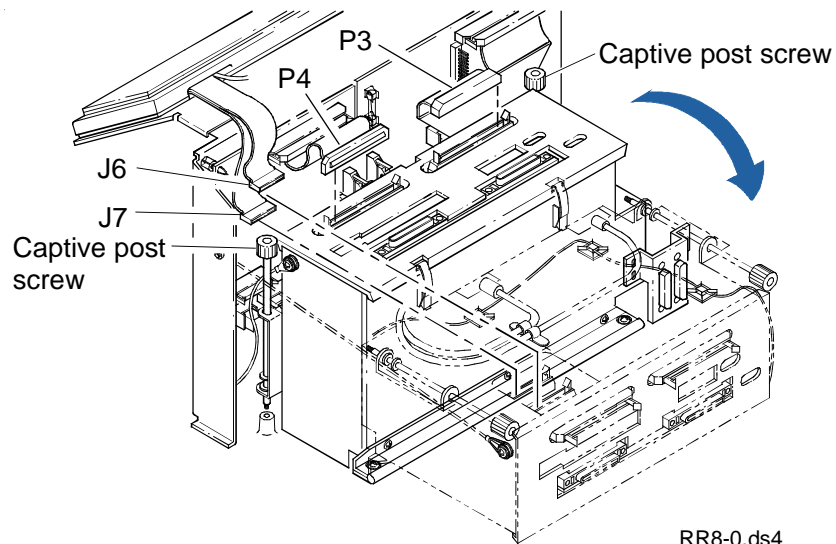

### 4 Removal & Replacement Procedures LCx<sup>®</sup> Analyzer Service Manual (Text P/N 47278-102)

### RR-8.0: Tilt Back Card Cage

4 - 88

PL-8

Chapter  
Contents

Manual  
Contents

CD Main  
Screen

Index