



VP-31: MATRIX COUNT CHECK (pg 1 of 1)


Time: 2 minutes

Purpose: To check cell counting ability

1. Load MEIA carousel ([Figure VP31](#)).

Blank Standard

Reaction
Cells

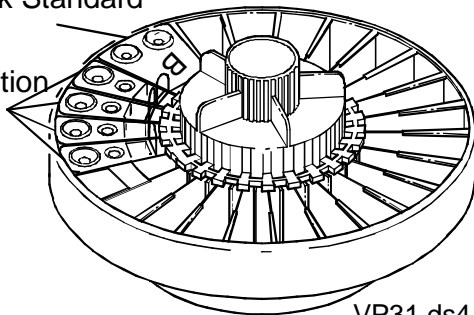

VP31.ds4

2. Begin Buffer Run:

- UTILITY
- OTHER
- OTHER
- SERVICE
- OTHER
- OTHER
- DIL_RUN

Assay # (assay the customer is using)

- STORE

3. Watch display.

4. When disposable count is displayed, press:

- STOP