



RR-5.7: LIQUID LEVEL SENSE BOARD (pg 1 of 2)

Removal



1. Remove Enclosure Assembly ([RR-1.3](#)).
2. Remove Liquid Level Sense Board:
 - a. Remove 1 mounting screw.
 - b. Slide board out.
 - c. Disconnect Liquid Level Sense Cables P1 thru P5.

Replacement


1. Install Liquid Level Sense Board:
 - a. Connect LLS cables P1 thru P5.
 - b. Slide board into guides on Power Supply.
 - c. Secure board to Power Supply with mounting screw.
2. Install Enclosure Assembly.

Verification

Perform these procedures in this order:

1. Level Sense Adjustment ([VP-12](#))
2. Dispense Check ([VP-13](#))

RR-5.7: Liquid Level Sense Board (pg 2 of 2)

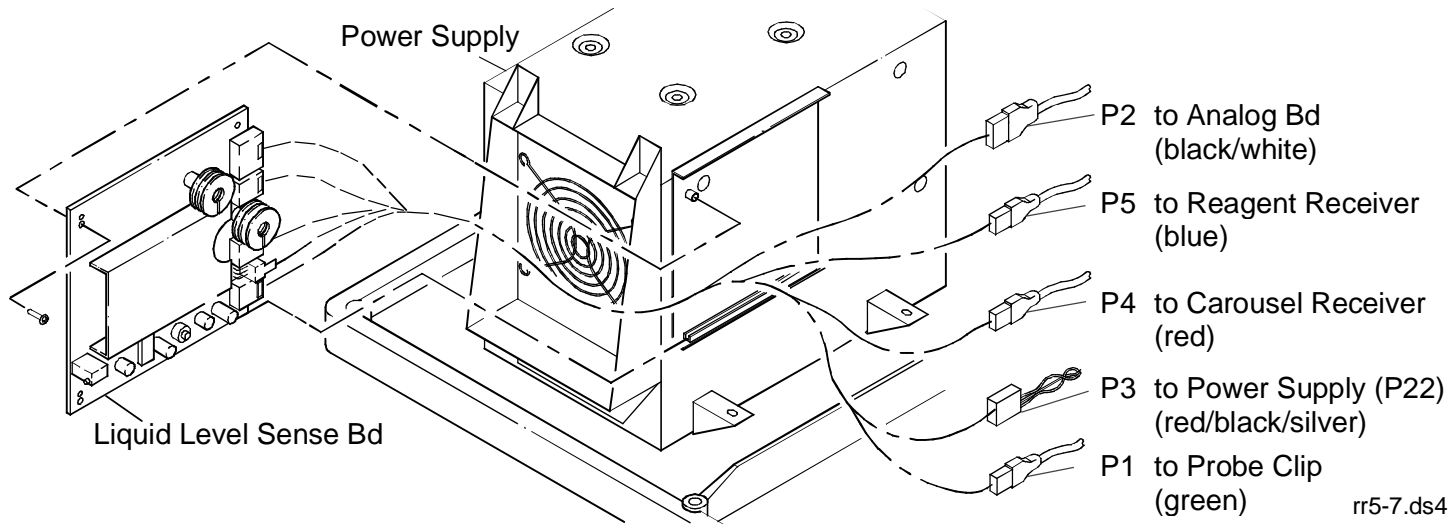

4 Removal & Replacement Procedures LCx[®] Analyzer Service Manual (Text P/N 47278-102)

RR-5.7: Liquid Level Sense Board

4 - 64

PL-5

Illustr. of Liquid
Level Sense Bd.

Chapter
Contents

Manual
Contents

CD Main
Screen

Index