


## VP-32: MATRIX DETECTION (pg 1 of 2)



Time: 2 minutes

Purpose: To verify that Matrix Detect Voltage is at the proper value

1. Load MEIA Carousel (Figure VP32).

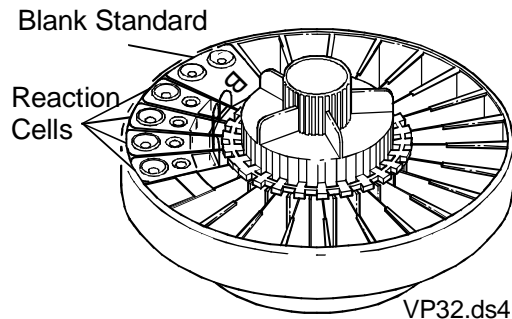

2. Close the door.
3. Print out System File 31:
  - SYSTEM
  - FILES
  - 31
  - PRINT
4. Enter MEIA Optics Hand Controls and adjust High Voltage:
  - MONITOR
  - HND\_CNTL
  - OPTICS
  - SET HVEnter MTRX DET VOLT 31.29 from printout.
  - STORE
  - LAMP (Verify that lamp is ON.)

## VP-32: Matrix Detection (pg 2 of 2)

5. Close door.
6. Wait for P= value to stabilize.
7. Press:
  - INDEXRecord the P= values for positions 1-8.
8. Verify specifications for the P= values.
  - Check 1: P= values for positions 1-4  
> parameter 31.30 value.
  - Check 2: P= values for positions 5-8  
< parameter 31.30 value.
  - Check 3: Lower value of positions 1 & 2  
minus the higher value of  
positions 7 & 8  $\geq 160$ .