

Chapter 4

Removal and Replacement Procedures

CHAPTER TABLE OF CONTENTS

RR Procedures (numeric order list)	4- 2
RR Procedures (alphabetical order list)	4- 4
Chapter Overview	4- 7
Removal & Replacement procedures	4- 9

WARNING !

Procedures in this chapter require following biohazard, electrical hazard, and electrostatic discharge precautions.



WARNING !

Tools must be decontaminated before and after servicing the analyzer. Tools must be decontaminated before they are returned to the tool kit.

CAUTION !

You must perform [VP-60: Decontamination](#) before you begin and after you complete any procedure that requires removing the analyzer cover.

REMOVAL & REPLACEMENT PROCEDURES (numeric order) (pg 1 of 2)

RR-1.3	Enclosure Assembly	4- 9	RR-3.9	Boom Assembly	4- 33
RR-2.3	Reagent Shroud	4- 12	RR-3.13	Air Thermistor Assembly	4- 36
RR-2.7	Protective Wire Cover (Interior Cover)	4- 13	RR-3.16	Air Heater Assembly	4- 38
RR-2.8	Air Duct Cover	4- 14	RR-3.18	Carousel Receiver Assy	4- 41
RR-2.11	Amplification Vial Retainer	4- 16	RR-3.19	Retainer Clip Assembly	4- 43
RR-3.1	Fan & Gasket Assembly	4- 19	RR-4.1	MEIA Optics Assembly	4- 45
RR-3.2	Reagent Receiver/ Heater Block	4- 22	RR-4.4	Lamp Power Supply	4- 48
RR-3.3	Carousel Motor Assembly	4- 25	RR-4.5	PMT HV Power Supply	4- 50
RR-3.4	Liquid Heater/Bracket Assy	4- 28	RR-4.6	Carousel Lock Sensor Cable Assembly	4- 52
RR-3.5	Probe Assembly	4- 31	RR-4.11	MEIA Lamp	4- 54
			RR-4.14	Remote Air Thermistor Assy	4- 57

Chapter Contents	Manual Contents	CD Main Screen	Index
---------------------	--------------------	-------------------	-------

RR Procedures (numeric list) (pg 2 of 2)

RR-5.1	Pump Assembly	4- 59
RR-5.3	Multivalve Block	4- 61
RR-5.7	Liquid Level Sense Board	4- 63
RR-5.8	Buffer Sensor Cable	4- 65
RR-5.10/12	Syringes	4- 67
RR-6.3	Printer Assembly	4- 70
RR-6.13	Main Power Supply	4- 72
RR-7.1/10	CPU/Digital I/O Boards	4- 75
RR-7.4	RS-232 Cable/Bracket Assembly	4- 79
RR-7.7	Card Cage Assembly	4- 81
RR-7.9	NOVRAM Replacement	4- 84

RR-8.0	Tilt Back Card Cage	4- 87
RR-8.1	Display Assembly (Touchpad)	4- 89
RR-8.1.1	Raise Display Procedure	4- 91
RR-8.3	Analog Board	4- 92
RR-8.5	Motor Driver Board	4- 94

Chapter Contents	Manual Contents	CD Main Screen	Index
------------------	-----------------	----------------	-------

REMOVAL & REPLACEMENT PROCEDURES (alphabetical order) (pg 1 of 3)

Air Duct Cover	RR-2.8	4- 14	Cable (Buffer Sensor)	RR-5.8	4- 65
Air Heater Assembly	RR-3.16	4- 38	Cable (Carousel Lock		
Air Thermistor Assembly	RR-3.13	4- 36	Sensor) Assembly	RR-4.6	4- 52
Air Thermistor (Remote)			Card Cage Assembly	RR-7.7	4- 81
Assembly	RR-4.14	4- 57	Card Cage (Tilt Back)	RR-8.0	4- 87
Amplification Vial Retainer	RR-2.11	4- 16	Carousel Lock Sensor		
Analog Board	RR-8.3	4- 92	Cable Assembly	RR-4.6	4- 52
Board (Analog)	RR-8.3	4- 92	Carousel Motor Assembly	RR-3.3	4- 25
Board (CPU/Digital I/O)	RR-7.1/7.10	4- 75	Carousel Receiver		
Board (Digital I/O)	RR-7.10	4- 75	Assembly	RR-3.18	4- 41
Board (Motor Driver)	RR-8.5	4- 94	Cover (Air Duct)	RR-2.8	4- 14
Boom Assembly	RR-3.9	4- 33	Cover (Protective Wire)	RR-2.7	4- 13
Buffer Sensor Cable	RR-5.8	4- 65	CPU/Digital I/O Boards	RR-7.1/7.10	4- 75

Chapter Contents	Manual Contents	CD Main Screen	Index
---------------------	--------------------	-------------------	-------

RR Procedures (alpha list) (pg 2 of 3)

Digital/I/O Board	RR-7.10	4- 75	Lamp Power Supply	RR-4.4	4- 48
Digital I/O/CPU Boards	RR-7.1/7.10	4- 75	Liquid Heater and		
Display Assembly			Bracket Assembly	RR-3.4	4- 28
(Touchpad)	RR-8.1	4- 81	Liquid Level Sense Board	RR-5.7	4- 63
Display (Raise Procedure)	RR-8.1.1	4- 91	Main Power Supply	RR-6.13	4- 72
Enclosure Assembly	RR-1.3	4- 9	MEIA Lamp	RR-4.11	4- 54
Fan & Gasket Assembly	RR-3.1	4- 19	MEIA Optics Assembly	RR-4.1	4- 45
Heater (Air) Assembly	RR-3.16	4- 38	Motor Driver Board	RR-8.5	4- 94
Heater Block/			Multivalve Block	RR-5.3	4- 61
Reagent Receiver	RR-3.2	4- 22	NOVRAM Replacement	RR-7.9	4- 84
Heater (Liquid) and			PMT HV Power Supply	RR-4.5	4- 50
Bracket Assembly	RR-3.4	4- 28	Power Supply (Lamp)	RR-4.4	4- 48

Chapter Contents	Manual Contents	CD Main Screen	Index
---------------------	--------------------	-------------------	-------

RR Procedures (alpha list) (pg 3 of 3)

Power Supply (Main)	RR-6.13	4- 72	Retainer Clip Assembly	RR-3.19	4- 43
Power Supply (PMT HV)	RR-4.5	4- 50	Remote Air Thermistor		
Printer Assembly	RR-6.3	4- 70	Assembly	RR-4.14	4- 57
Probe Assembly	RR-3.5	4- 31	RS-232 Cable/		
Protective Wire Cover	RR-2.7	4- 13	Bracket Assembly	RR-7.4	4- 79
Pump Assembly	RR-5.1	4- 59	Shroud (Reagent)	RR-2.3	4- 12
Raise Display Procedure	RR-8.1.1	4- 91	Syringes	RR-5.10/5.12	4- 67
Reagent Receiver/ Heater Block	RR-3.2	4- 22	Tilt Back Card Cage	RR-8.0	4- 87
Reagent Shroud	RR-2.3	4- 12	Touchpad Assembly (Display)	RR-8.1	4- 89
Receiver (Carousel) Assembly	RR-3.18	4- 41			
Receiver (Reagent)/ Heater Block	RR-3.2	4- 22			

Chapter Contents	Manual Contents	CD Main Screen	Index
----------------------------------	---------------------------------	--------------------------------	-----------------------

CHAPTER OVERVIEW

This chapter provides:

- step-by-step procedures for removing and replacing LCx® Analyzer parts
- illustrations to accompany the procedures


Links

- RR procedure and its related illustration
- RR and related Parts List
- RR and the alphabetical and numerical lists of RR procedures
- Verification Procedures (VPs) in the RR and the related VPs
- Board and Card Cage illustrations

Numbering Scheme

The Removal/Replacement procedures in this chapter are indexed by number to the related Parts List in Chapter 3.

Example: RR-2.5 is the Removal/ Replacement procedure for Item 5 on Parts List PL-2.

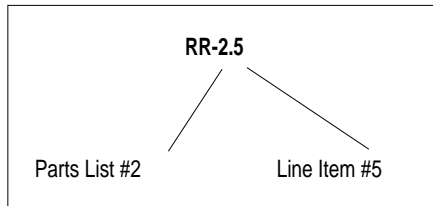

This page left intentionally blank.

Chapter Contents	Manual Contents	CD Main Screen	Index
----------------------------------	---------------------------------	--------------------------------	-----------------------