

RR-5.8: BUFFER SENSOR CABLE

(pg 1 of 2)



Removal

1. Remove Enclosure Assembly ([RR-1.3](#)).
2. Remove Buffer Tray and 2 mounting screws.
3. Remove Buffer Sensor Cable:
 - a. Remove Cable Retainer.
 - b. Remove Buffer Sensors.
 - c. Disconnect Buffer Sensor Cable connector.


Replacement

1. Install Buffer Sensor Cable:
 - a. Secure Buffer Sensors to Buffer Tray with 2 mounting screws.
 - b. Install Cable Retainer.
 - c. Connect Buffer Sensor Cable connector.
2. Install Buffer Tray. Secure tray over standoffs with 2 mounting screws.
3. Install Enclosure Assembly.

Verification

Perform Buffer Sensor Check ([VP-15](#)).

RR-5.8: Buffer Sensor Cable (pg 2 of 2)

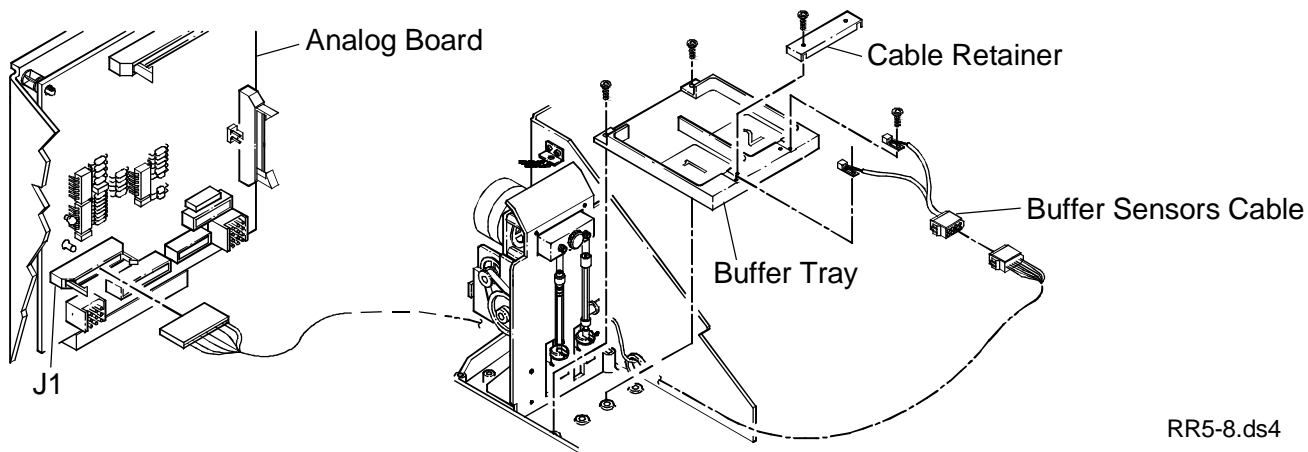

4 Removal & Replacement Procedures LCx[®] Analyzer Service Manual (Text P/N 47278-102)

RR-5.8: Buffer Sensor Cable

4 - 66

PL-5

Chapter
Contents

Manual
Contents

CD Main
Screen

Index