

 ABBOTT ADD	INDEX TECHNICAL SERVICE BULLETIN
---	---

PRODUCT: ALCYON (TM) 300 ISE (122)	DATE: 23-SEP-98
--	---------------------------

TSB #	IMPLEMENTATION	SUBJECT	EFFECTIVITY DATE
122-002	N - SN 1831 & below	ALCYON (TM) Rev. 0 Software Release Ver. 261095/0142	20-APR-98
122-001	N - See TSB	Modified Wash Pump with New O-ring	07-OCT-97

PENDING -	TSB index number has been reserved for a future TSB.
CANCELLED -	TSB index number is cancelled.
INCORPORATED -	TSB was incorporated into another document or manual.
OBSOLETE -	TSB no longer applies.
COMPLETE -	TSB implementation is complete.

END OF DOCUMENT

 ABBOTT ADD	<h1>TECHNICAL SERVICE BULLETIN</h1>
---	-------------------------------------

SUBJECT:

ALCYON (TM) Rev. 0 Software Release Ver. 261095/0142TSB#: **122-002**ORIGINATOR: **Emile Diou/Bernard Dabosi**APPROVED: **Christie McCain 05/21/98**

PRODUCT:

ALCYON 300 ISE (122)REF. ECN: **Alcyon France # I-11**Trademark: **ALCYON™** is a trademark of Abbott Laboratories.

IMPLEMENTATION: <input checked="" type="radio"/> Immediate <input type="radio"/> Next Service Call <input type="radio"/> Next Failure <input type="radio"/> Optional Instruments Requiring Modification: SN 1831 and below, Abbott distributed instruments only	TSB Part/Kit #: <u>A201300A</u> TSB Effectivity/ Part(s) Availability: <u>20-APR-98</u> TSB Tracking by Serial # required (IMMEDIATE TSB's ONLY) <input checked="" type="radio"/> YES <input type="radio"/> NO	Upgrade Time: <u>0.5 Hr.</u> Validation Time: <u>0.5 Hr.</u> Total Mod. Time: <u>1.0 Hr.</u> **NOTE** The instrument must be at TSB Level <u>n/a</u> prior to performing this TSB.
---	--	---

I. DISTRIBUTION:

Worldwide

****Potential Biohazard & Voltage Hazard. Observe Proper Safety Precautions.****

II. PURPOSE:

The purpose of this TSB is to inform the field of the new software released to resolve the issue of a possible mismatch of test results with the Sample IDs.

Problem: When four or more samples are being processed, if samples selected for Run from the worklist are not selected in the sequence in which they appear in the worklist, some test results can be reported under an incorrect Sample ID.

Solution: The new software has been designed to resolve the problem. There are no changes to the user interface, only to the software code.

III. ADMINISTRATIVE NOTES:

NOTE: **Since this TSB is IMMEDIATE, tracking is required.**

International: This software package will be distributed to ACS in Delkenheim for distribution to Europe, Africa and the Middle East; to Singapore for distribution to SMITV; and to each Country Service Manager outside those areas. Each Service Manager will decide whether to send this package directly to the Customer or have an Abbott representative perform the installation. When the installation of the new software is complete, FAX the "Software Installation Report" printout to:

Alcyon Customer Systems Engineering
Dallas, Texas USA
Attn: Emile A. Diou

****Potential Biohazard & Voltage Hazard. Observe Proper Safety Precautions.****

Fax Number: +1-972-753-3525

Additionally....

For instruments in Europe, Middle East, and Africa, please copy

Amarjit Panesar @ +49-6122-58-1521

For instruments in Southeast Asia and the Far East, please copy

Heng Chiang Wee @ +65-270-8873

For instruments in Latin America, please copy

Mayra Vazquez @ 1+847-938-6027

Mark off the #2 on the TSB Modification Control Sticker after the analyzer has been upgraded.

IV. SPECIAL TOOLS:

None.

V. PARTS:

1 ea. - Abbott ALCYON Main Program Diskette, Ver. 261095/0142: P/N A201300A, Abbott LN 05D64-01

Sufficient quantities of the following will be sent to each country where instruments have been shipped:

- 1) The new Main Program Diskette
- 2) Cover letter for the Country Service Managers
- 3) Official ADD Customer Letter (explaining change to Customer)
- 4) Installation/Upgrade Instructions

****Potential Biohazard & Voltage Hazard. Observe Proper Safety Precautions.****

REPLACED PARTS:

Discard any older versions of the ALCYON Main Program software.

COMPATIBILITY:

Compatible with all ALCYON Analyzers produced for Abbott.

VI. PROCEDURE:

NOTE: Make sure the printer is on and connected before performing this procedure.

**NOTE: If any problems occur during this process, DO NOT try to install the software again.
Contact your Area Specialist or WW Customer Systems Engineer.**

1. Print the current Configuration parameters from the current software version.
2. Select CONFIGURATION/TOOLBOX from the main menu.
3. Select PRINT OPTIONS, then press SHIFT key + PRINT SCREEN key to print.
4. Select OTHER CONFIGURATION PARAMETERS, then press SHIFT key + PRINT SCREEN key to print.
5. Select INTERFACE, then INTERFACE SETUP, then press SHIFT key + PRINT SCREEN key to print.
6. Select SAMPLE CODEBAR READER SETUP, then press SHIFT key + PRINT SCREEN key to print.

****Potential Biohazard & Voltage Hazard. Observe Proper Safety Precautions.****

7. Switch the instrument OFF. Put the Saving Methods Diskette into the floppy disk drive. Switch the instrument ON and perform a backup of the Methods Savings following the procedure in the Operations Manual. When complete, remove the diskette and switch the instrument OFF.
8. Insert the new system software diskette into the disk drive, and switch the instrument ON. The system will boot-up to the Installation Language Menu.
9. Select a language and press ENTER. Follow the screen instructions by completing the fields (Installer, Laboratory, and Town are required fields). Press F10 when complete. The message, "Printing Report you must Send" will appear and the system will print an installation report. The system will then display "Copying Files" while the software is loading.

Please FAX the printed report to the person(s) indicated in the Administrative Notes above.

10. Remove the system software diskette when prompted on the display. Then, select CONFIRM, and press ENTER. The instrument will initialize.
11. Verify the Configuration information by comparing the printed information from **Steps 2 through 6**. Make changes if necessary. Select F10 to validate any changes.
12. Verify the Customer's Methods Definitions to be sure they have not changed.
13. Store the new System Software diskette in a safe location for future use. Discard any diskettes containing older versions of the ALCYON Main Program Software.

MODIFICATION STEPS:

Refer to PROCEDURE section above.

CHECKOUT:

Recycle power to the system then run assay (any) controls to assure the instrument is functioning properly.

MODIFICATION CONTROL STICKER UPDATE:

Mark off block #2 on the TSB Modification Control Sticker.

END OF DOCUMENT

 <p>ABBOTT ADD</p>	<h1>TECHNICAL SERVICE BULLETIN</h1>
--	---

SUBJECT:

Modified Wash Pump With New O-RingTSB#: **122-001**ORIGINATOR: **Gary Tompkins/Bernard Dabosi**

PRODUCT:

APPROVED: **Christie McCain 9/21/98****ALCYON (TM) 300 ISE (122)**REF. ECO: **Alcyon, France M160**Trademark: **ALCYON (TM)** is a trademark of Abbott Laboratories

<p>IMPLEMENTATION:</p> <div style="border: 1px solid black; padding: 5px; background-color: #f0e68c;"> <input type="radio"/> Immediate <input checked="" type="radio"/> Next Service Call <input type="radio"/> Next Failure <input type="radio"/> Optional </div> <p>Instruments Requiring Modification: S/Ns 1613, 1614, and 1639 through 1646</p>	<p>TSB Part/Kit #: <u>A072600B</u></p> <p>TSB Effectivity/ Part(s) Availability: <u>07-OCT-97</u></p> <p>TSB Tracking by Serial # required (IMMEDIATE TSB's ONLY)</p> <div style="border: 1px solid black; padding: 5px; background-color: #f0e68c;"> <input type="radio"/> YES <input checked="" type="radio"/> NO </div>	<p>Upgrade Time: <u>.5 Hrs.</u></p> <p>Validation Time: <u>.5 Hrs.</u></p> <p>Total Mod. Time: <u>1.0 Hrs.</u></p> <p>**NOTE** The instrument must be at TSB Level <u>n/a</u> prior to performing this TSB.</p>
--	--	---

****Potential Biohazard & Voltage Hazard. Observe Proper Safety Precautions.****

I. DISTRIBUTION:

Worldwide except United States.

II. PURPOSE:

Assembly methods caused premature failure of the exit port fittings of the Wash Pump. This TSB directs the replacement of the entire Wash Pump, which has added o-rings to eliminate leaks and broken fittings.

III. ADMINISTRATIVE NOTES:

This TSB is to be incorporated as a NEXT SERVICE CALL upgrade of the Wash Pump, Alcyon Ref. Number A072600B (new number).

This modification should be performed in the course of normal service.

Europe: Area Customer Service, Delkenheim, should communicate to WWCS in Dallas, FAX# 972-753-3525, when the instruments in Europe and surrounding areas have been upgraded.

ROW: Area Service Managers should communicate to WWCS in Dallas, FAX# 972-753-3525, when the instruments in their respective areas have been upgraded.

IV. SPECIAL TOOLS:

No special tools required. (Metric tools are considered standard for the ALCYON Analyzer).

V. PARTS:

REPLACED PARTS:

The pumps replaced during the modification should be scrapped according to country policy.

COMPATIBILITY:

This modification is compatible with all existing Rev. 0 ALCYON Analyzers.

Any spare Wash Pump ordered as of the effectivity date of this TSB will be modified with the new o-ring and assembly techniques.

VI. PROCEDURE:

MODIFICATION STEPS:

Set-up: **Be sure to wear proper safety equipment (gloves, glasses and lab coat).**

1. Remove the water input lines & sensor from the water supply tank and short the sensors together. Perform a probe wash for 1 minute to air purge the pump and wash lines. Select
 - System Configuration/Tool Box
 - Washing Probes
 - ENTER to Begin and End
2. Power the system OFF using the key switch, then the main power switch and unplug the power cord. Place paper towels appropriately under the pump to catch any liquid when the tubing is disconnected.

****Potential Biohazard & Voltage Hazard. Observe Proper Safety Precautions.****

3. Disconnect the inlet tubing from the top of each pump head. The barbed fittings must remain, they are secured by glue. Remove the Check Valve and save it.
4. Unscrew the output Teflon tubing from the bottom of each pump head. The white fittings must remain, they are secured by glue.
5. Disconnect the pump power cable at connector M2.
6. Using a 8mm hex socket wrench, loosen the four screws which secure the pump and remove the pump. Retighten the four shock mounts as necessary.
7. Install the new pump by following Steps 3 through 6 in reverse order. Be sure to correctly attach the Check Valve and output and input lines.

VALIDATION:

Turn the instrument ON and enter into the Reglages software. Be sure the supply tubes are in the water supply tank and the tank has sufficient distilled water in it.

1. At the Reglages main menu, select
 - J, HYDRAULIC submenu, then select
 - A, PUMP ON


Allow to run for 1 minute to fill the hydraulic circuit while verifying there are no leaks at the pump fittings and tubing connections. Select

- B, PUMP OFF

2. Pump Flow Check
 - a) Remove the Reagent Probe from the Pre-heater but DO NOT disconnect the tubing from the probe. (Be sure to loosen thumb screw enough so the probe is not scratched while removing.)
 - b) Place the probe in a clean, empty container with a capacity of no less than 50ml (>= 50ml).

****Potential Biohazard & Voltage Hazard. Observe Proper Safety Precautions.****

- c) Remove only the one large waste tubing (A072300A) from the waste bottle and place it in an empty container with capacity of greater than 100ml (>100ml). This tube will be used to check the output of the sample head (combined from Sample Probe and wash well lines).

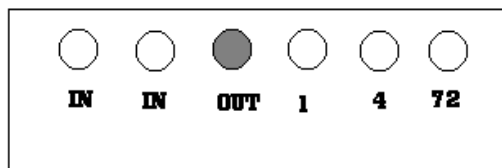

Waste tubing comes from this port (OUT) at rear of the ALCYON 300 ISE instrument

- d) At the Reglages main menu, select
- J, HYDRAULIC submenu, then select
 - A, PUMP ON to run the pump FOR EXACTLY 15 SECONDS. Then select
 - B, PUMP OFF.
- e) Measure the volume in the 'reagent' container. The target volume is **25ml**; the range is 25ml to 30 ml.
- f) Measure the volume in the 'waste' container. The target volume is **60ml**; the range is 60ml to 80ml.
- g) Return the waste tubing to waste container and the reagent probe to the probe holder when the check is complete and acceptable.

If measured volume is not correct, perform the Pump Flow Adjustment, Step 3.

3. Pump Flow Adjustment

This procedure must be performed if the volumes measured during the Flow Check were not

****Potential Biohazard & Voltage Hazard. Observe Proper Safety Precautions.****

correct.

a) If incorrect volume from reagent head:

- If the volume from the reagent head is less than 25ml, turn the reagent flow adjustment screw (see **Figure 1**) in the counter-clockwise direction for 1/16 of a turn to increase the flow. Repeat the Flow Check.
- Perform the Solenoid Valve Check to verify that the valve does not leak. If the volume from the reagent head is greater than 25ml and less than 30ml but the solenoid valve leaks, turn the reagent flow adjustment screw (see **Figure 1**) in the clockwise direction for 1/16 of a turn to decrease the flow to 25ml. The pressure may be too high, causing the valve to leak. (More than one adjustment may be necessary.) Repeat the Flow Check.
- Repeat the Solenoid Valve Check. If the reagent solenoid valve continues to leak and the volume is right at the target value of 25ml, replace the valve (A143600B) then repeat all the checks.

b) If incorrect volume on pump sample head:

- If the volume from the sample head is less than 60ml, turn the sample flow adjustment screw (see **Figure 1**) in the counter-clockwise direction for 1/16 of a turn to increase the flow. Repeat the Flow Check.
- Perform the Solenoid Valve Check to verify that the valve does not leak. If the volume from the sample head is greater than 60ml and less than 80ml but the solenoid valve leaks, turn the sample flow adjustment screw (see **Figure 1**) in the clockwise direction for 1/16 of a turn to decrease the flow. The pressure may be too high, causing the valve to leak. (More than one adjustment may be necessary.) Repeat the Flow Check then the Solenoid Valve Check.
- Repeat the Solenoid Valve Check. If the sample solenoid valve continues to leak and the volume is right at the target value of 60ml, replace the valve (A143500B) then

****Potential Biohazard & Voltage Hazard. Observe Proper Safety Precautions.****

repeat all the checks.

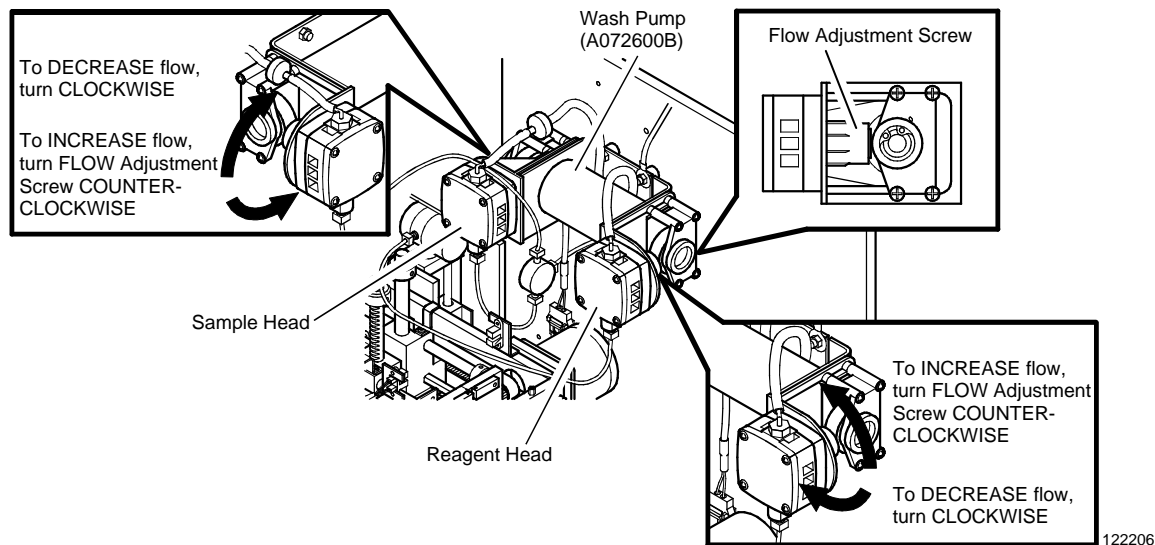


Figure 1

4. Solenoid Valve Check

This step should be repeated after any adjustment is made to the pump flow.

a) Reinstall the Reagent Probe into the Pre-heater.

****Potential Biohazard & Voltage Hazard. Observe Proper Safety Precautions.****

b) Manually move the Reagent and Sample Transfer Arm to the UP positions.
Keep them above the wash well.

c) At the Reglages main menu, select

- J, HYDRAULIC submenu, then select
- A, PUMP ON

Allow the pump to run for 3 minutes. Check for a proper probe stream in both probes. Continue to observe the pump fittings and tubing for leaks.

d) After 3 minutes, alternate opening and closing the valves every 2 seconds for a period of 10 seconds as follows:

- D, ELECTROVALVES CLOSING for 2 seconds
- C, ELECTROVALVES OPENING for 2 seconds

Repeat this 5-10 times. The probes MUST NOT LEAK when the pump is running and the VALVES are CLOSED.

[If a probe still leaks, check the appropriate circle Ts, syringe, and tubing. If probe continues to leak, replace the corresponding solenoid valve.]

e) Select

- PUMP OFF to turn the pump off.

5. Exit the Reglages software

- ESC
- ESC
- Y (YES)
- ENTER key to exit to the C:\> prompt

6. Restart the system using the key switch.

7. Complete the current service call.

****Potential Biohazard & Voltage Hazard. Observe Proper Safety Precautions.****

CHECKOUT:

See VALIDATION above.

MODIFICATION CONTROL STICKER UPDATE:

Mark off the #1 on the modification control sticker to indicate that TSB 122-001 has been incorporated.

END OF DOCUMENT