
## VP-19: CAROUSEL MOTOR AND GEAR ADJUSTMENT (pg 1 of 2)

Time: 10 minutes

## Purposes:

- To set carousel gear position on carousel motor shaft
- To position motor at correct height to ensure positive carousel drive
- 1. Adjust carousel gear positioning (Figure VP19):
  - a. Place Carousel Motor Gear Alignment Tool in analyzer with cutout over gear.
  - b. Loosen set screw.
  - c. Center gear in cutout.
  - d. Tighten set screw.

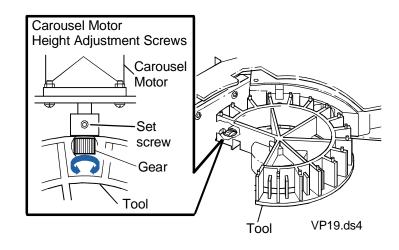

5 Verification Procedures VP-19: Carousel Motor and Gear Adjustment 5 - 50 LCx® Analyzer Service Manual (Text P/N 47278-102)

Chapter	Manual	CD Main	Index
Contents	Contents	Screen	

## VP-19: Carousel Motor and Gear Adjustment (pg 2 of 2)

- 2. Adjust carousel motor height:
  - a. Place Tool in analyzer with cutout over gear.
  - b. Put downward pressure on Tool and rotate Tool clockwise (CW) and counterclockwise (CCW) over gear. Teeth must fully mesh, but Tool must not be raised by gear.
  - c. Loosen Carousel Motor Height Adjustment Screws.
  - d. Adjust motor up or down until Tool rotates properly over gear.

5	Verification Procedures	VP-19: Carousel Motor and Gear Adjustmen	nt 5 - 52
	LCx® Analyzer Service Manual (Text P/N 47278-	02)	

Chapter	Manual	CD Main	Index
Contents	Contents	Screen	