

RR-8.3: ANALOG BOARD (pg 1 of 2)



Removal

1. Remove Enclosure Assembly ([RR-1.3](#)).
2. Tilt back Card Cage ([RR-8.0](#)).
3. Disconnect 9 Analog Board cables.
4. Remove Analog Board.

Replacement

1. Place Analog Board into board guide; press into place.
2. Connect 9 Analog Board cables.
3. Reposition and secure Card Cage.
4. Install Enclosure Assembly.

Verification


Perform these procedures:

- Dispense Check ([VP-13](#))
- MEIA Temperature Cal ([VP-27](#))
- MEIA Photo Check ([VP-26](#))

PL-8	Block diagram: Analog Bd.
------	------------------------------

Chapter Contents	Manual Contents	CD Main Screen	Index
---------------------	--------------------	-------------------	-------

RR-8.3: Analog Board (pg 2 of 2)

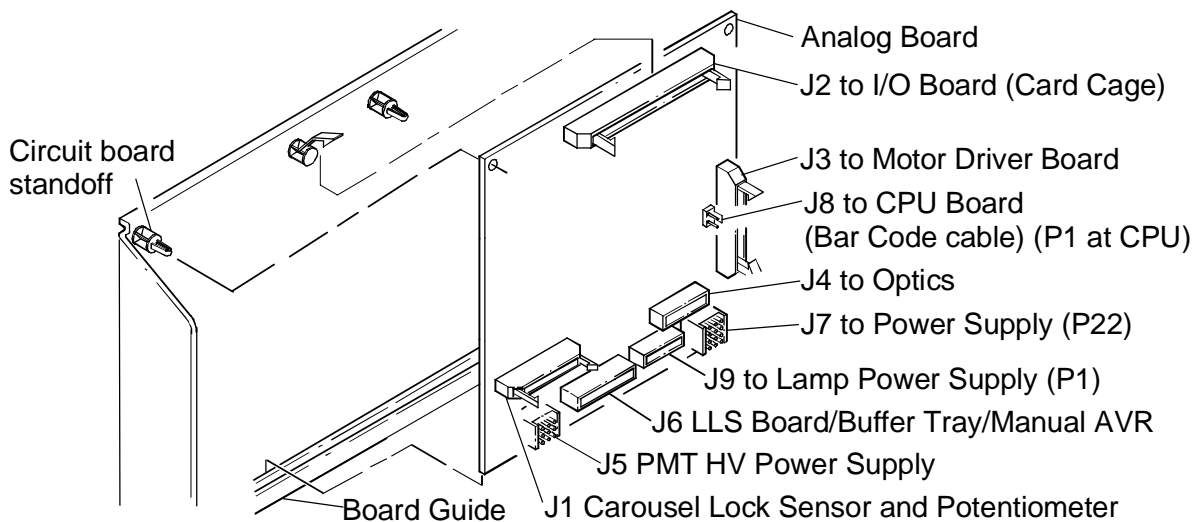

RR8-3.ds4

4 Removal & Replacement Procedures LCx[®] Analyzer Service Manual (Text P/N 47278-102)

RR-8.3: Analog Board

4 - 93

PL-8

Block diagram:
Analog Bd.

Chapter
Contents

Manual
Contents

CD Main
Screen

Index